


# Nucleic Acids / Disease Biomarker Profiling II

## Conditions


|               |                                                            |
|---------------|------------------------------------------------------------|
| Column:       | ACE 3 C18-PFP<br>ACE 3 C18                                 |
| Dimensions:   | 100 x 4.6 mm                                               |
| Part Number:  | ACE-1110-1046, ACE-111-1046                                |
| Mobile Phase: | 33 mM potassium phosphate pH 6.2 with KOH/MeOH (88:12 v/v) |
| Flow Rate:    | 1 mL/min                                                   |
| Temperature:  | Ambient                                                    |
| Detection:    | UV, 260 nm                                                 |

### Deoxyribonucleosides

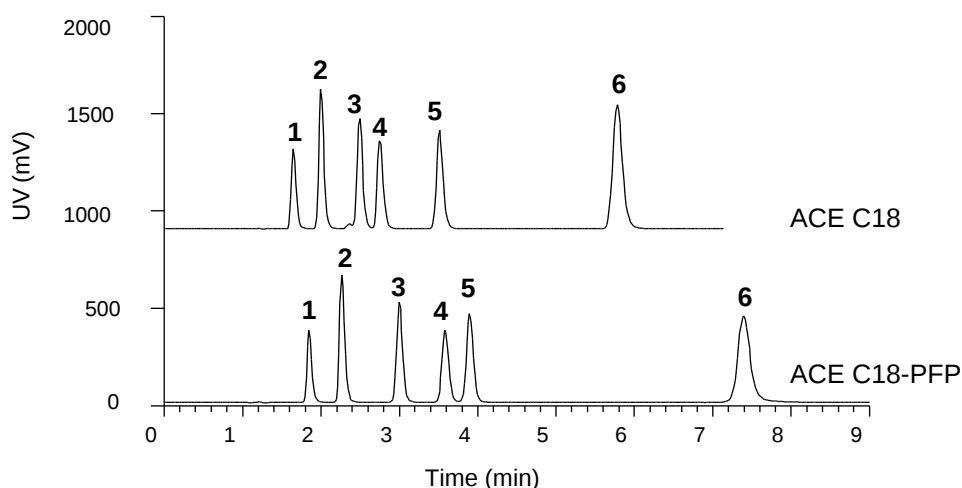

1. Deoxycytidine
2. Deoxyuridine
3. Deoxyinosine
4. Deoxyguanosine
5. Thymidine
6. Deoxyadenosine

### Ribonucleosides

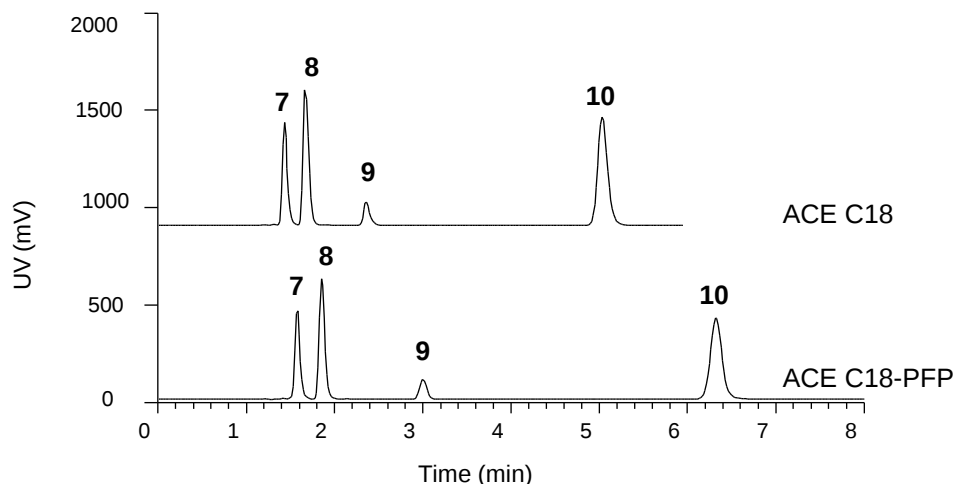

Reproduced with permission of Department of Medical Biochemistry and Biophysics, Umeå University, Sweden

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit  
[www.ace-hplc.com](http://www.ace-hplc.com) or email: [info@ace-hplc.com](mailto:info@ace-hplc.com)