

Sennosides in Traditional Chinese Medicine



Application #AN1390

Conditions


Column: ACE 3 C18-PFP
Dimensions: 150 x 4.6 mm
Part Number: ACE-1110-1546
Mobile Phase: A: 0.75% acetic acid in H₂O
B: MeCN/MeOH (90:10 v/v)

Time (mins)	%B
0	9
23	28
40	28

Flow Rate: 0.6 mL/min
Temperature: 35 °C
Detection: UV, 271 nm
Sample: Herbal tea bag containing Folium Sennae, Peppermint, Folium Mori, Folium Nelumbinis, Glycyrrhiza Uralensis and Lalang Grass Rhizome

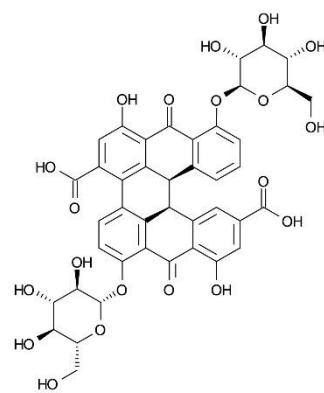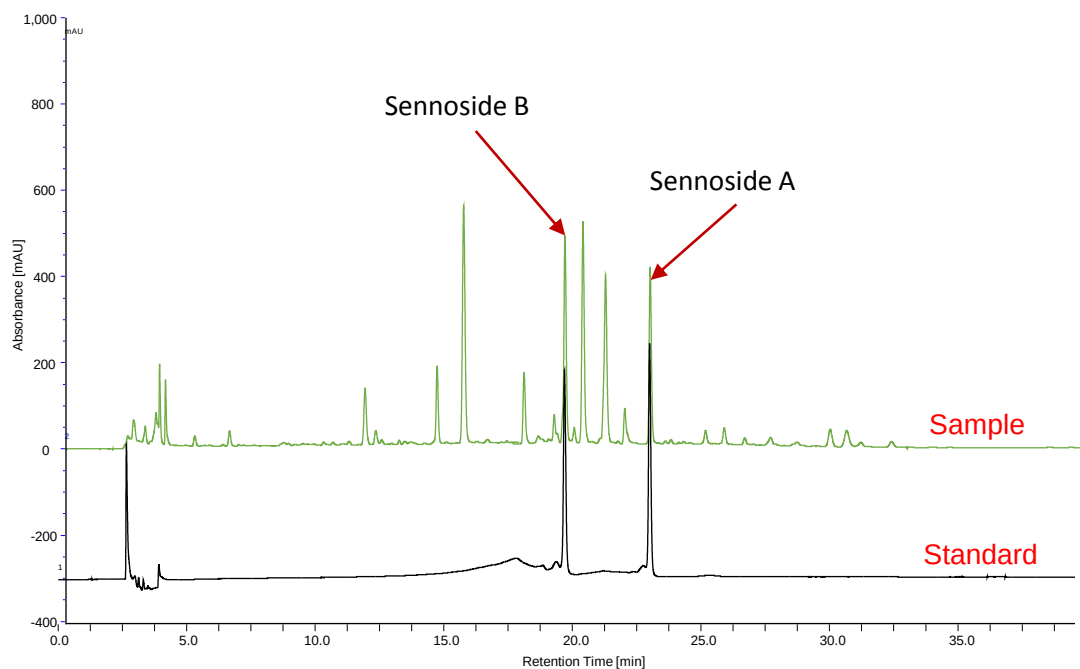

Sennoside A

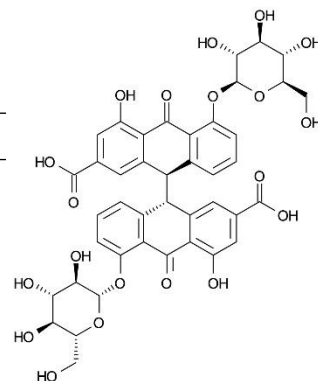

Sennoside B

Reproduced with permission of Vercotech Inc., Taiwan

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com